


THANK YOU!




REGENERATIVE SKIN THERAPY

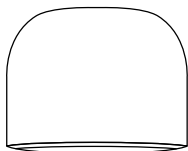


**Hygiene Cap**  
Protects the needles during transportation and reduces the entry of germs

**1.0 mm**  
Needle Length

**20 Needles**  
Hypoallergenic  
Sterile 24K  
Gold-Plated  
Titanium  
Surgical Needles

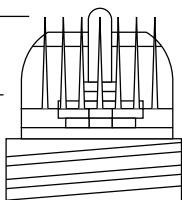

**5 ml**  
Vial Capacity

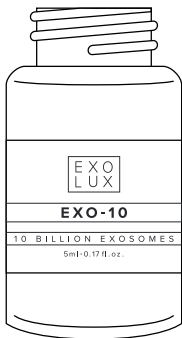

Do Not Use If This Package Has Been Opened Or Damaged.



KEEP OUT OF CHILD'S REACH



### Instructions:

1. Remove the applicator from the packaging.
2. Sanitize the top of applicator about 10 minutes with 75% Isopropyl Alcohol.
3. Wash your skin with warm water.
4. Put the essence of the liquid into the applicator vial. Carefully screw the top back onto the vial. Gently shake the vial 2-3 times.
5. Remove the safety and hygiene cover.
6. Pat the skin surface with the applicator gently until the essence of the liquid is fully absorbed by the face.

### Cautions:

1. Do not share applicator with another person for hygiene.
2. Always keep the applicator clean.
3. Keep the applicator out of reach of children.



KEEP OUT OF CHILD'S REACH



[WWW.EXO-LUX.COM](http://WWW.EXO-LUX.COM)